## HF2002

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Dear Customer,

we are pleased that you have decided to purchase a product from our range. Please read these operating instructions carefully and carefully and keep them accessible for other users. Observe the application instructions and notes.

With kind regards  
*the Novidion team*

## Hair wash basin HF2002 Operating instructions

### 1. Intended use

The function of the Pulox hair wash basin is to effortlessly wash the hair of people in need of care or who are bedridden while sitting or lying down, without the person having to be moved. The product is suitable for use in the private sector (home care) and in medical facilities (hospital, clinic, retirement home).

The device can be used by persons aged 6 years and over and should only be operated by persons who are in full possession of their physical, sensory and mental faculties.

### 2. Contraindications

None

### 3. Safety instructions

Please read the safety instructions carefully before using the product. Operators do not require professional training, but should only use this product after they have fully understood the requirements in these operating instructions.

1. If there is visible damage to the basin or accessories, please do not use the product and contact the manufacturer or distributor.
2. Only use the product in accordance with the operating instructions and not for any other purposes not described in these operating instructions.
3. This product is not a toy. Store the product and packaging out of the reach of children.
4. All serious incidents related to the device shall be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.
5. This product can be used several times. Please follow the cleaning instructions after each use.

### 4. Product characteristics

#### Product parameters:

Product name: Hair wash basin

#### Performance characteristics:

A. Easy to use and suitable for non-

Model: HF2002  
 Classification: Class I medical  
 device

professionals  
 B. Handy, lightweight and mobile  
 C. Effective aid for cleaning immobile  
 patients

Materials contained:

- PP (Polypropylene)
- PVC (Polyvinyl chloride)


## 5. Accessories

In addition to the inflatable pool with drain hose, the following components are also included:

- 1x Shower bag with hand shower
- 1x Air pump


## 6. Product structure

The main product consists of an inflatable PVC basin with drain and lockable valve. The basin can be inflated by mouth or with the supplied air pump. The shower bag (10 litre capacity) with suspension and hand shower can be filled via a valve. The 750 mm long hose with shower head ensures freedom of movement and flexible use.

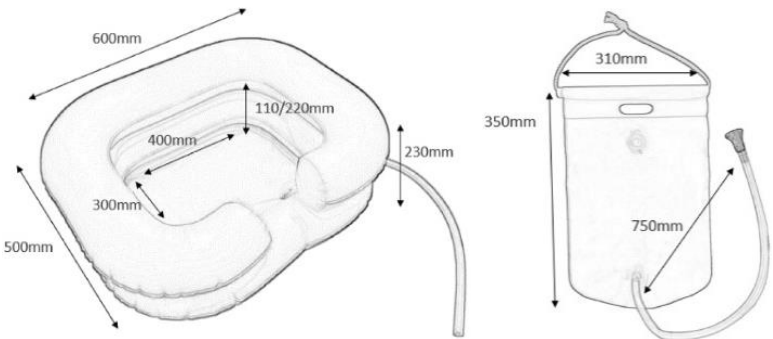
## 7. Application

1. Inflate the washbasin with the air pump and close the valve. Avoid over-inflating the product to prevent it from bursting. Do not under-inflate it, otherwise it could become unstable.
2. Place the washbasin in a flat and stable location.
3. Place the drain hose of the sink in a tank (e.g. a bucket) next to the sink so that the water can drain away.
4. Open the inlet valve of the shower bag, fill it as required and close the valve tightly. The bag has a reinforced handle and can be hung on the bed or similar fixtures near

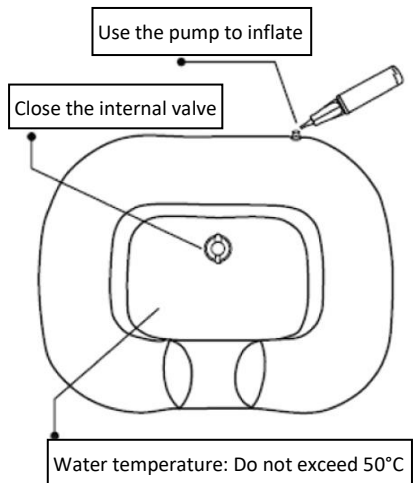
the place of use.

5. Place the patient with their head in the neck opening of the washbasin so that it is comfortable for them.
6. The long, handy hose has a shower head with an adjustable tap with which you can control the gentle water jet.
7. After use, open the valve on the inside of the washbasin so that the water can drain into the container provided.
8. Drain all the liquid and clean the basin and shower bag before allowing the components to dry completely.

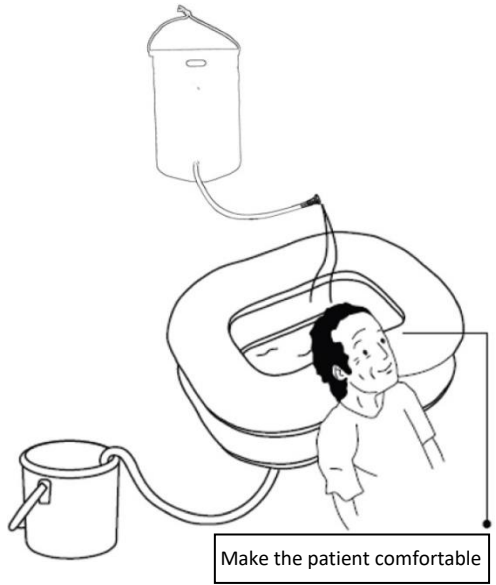Open the drain valve

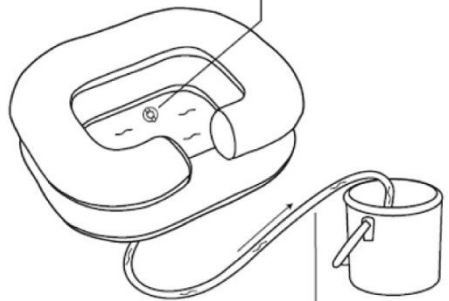Drain residual water into the tank

## 8. Storage conditions

Clean and dry all components carefully. Do not store the product in places where it is exposed to high temperatures or near sharp objects that could damage the material.

## 9. Applicable symbols

Symbol	Meaning	Symbol	Meaning
	Attention! Observe the accompanying document (user manual)		Attention! Observe the accompanying document (user manual)
	Medical device		Serial number
	Recyclable		Protect from direct sunlight
	Batch number		Store in a dry place
	Date of manufacture		Manufacturer
	European representative		Product complies with Directive 93/42/EEC of 14 June 1993 concerning medical devices

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Novidion GmbH
Foshan Hongfeng Co., Ltd
MedNet EC-REP GmbH